
**E-BODA** *life*

# Electronic body scale

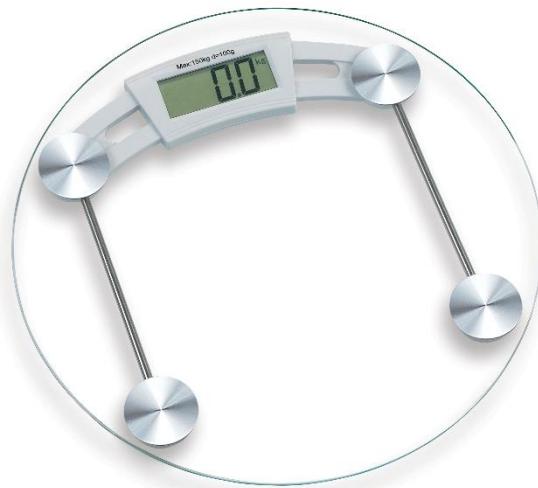**Carefully read the manual**

**CEP 1221**

### **Beware:**

- To avoid scratching the glass platform take off your street footwear or hard sole wear before weighting.
- The glass platform may be slippery if it is wet. Before weighting please make sure that the sole of your foot and the glass platform are dry.
- For cleaning you can use a fine towel lightly soaked in alcohol. Do not use water, avoid exposing the scale to extreme cold and heat.
- Do not jump on the scale, do not hit it or drop it. Its precision may deteriorate due to violent shocks.
- This scale is not built for professional use.
- We recommend replacing the battery if functioning errors occur. If the errors do not stop appearing, please call our customer service or schedule an appointment on [www.e-boda.ro](http://www.e-boda.ro).

### **Functionality:**

- High precision system
- Glass platform
- Capacity ranging from 6Kg to 150 Kg
- 100g divisions
- LCD screen
- Resetting to 0
- Automatically turns on
- Automatically turns off
- Low battery indicator
- Overload indicator
- Switch to change in 3 measuring units: kg, lb, st.
- Recommended to persons ranging from age 10 to 80
- Uses one 3V CR 2032 battery

**CEP 1221**

### **Mode of operation:**

- Mount the battery in the pack inside the assigned space on the back of the scale respecting the marked polarity on the battery and on the cover of battery slots.
- Place the scale on a rigit, perfectly plain and clean surface.
- Press on the glass platform and the scale will show a series of number "8" and then 0.0Kg.
- Gently step on the scale with both legs. Place your feet closer to the center of the glass platform, then straighten up and do not move for a few seconds. The screen will show your real weight after it clips three times.
- The screen will automatically turn off if the scale is not used anymore.

### **Caution:**

- For accurately measuring your weight we recommend weighting yourself a few days in a row at the exact same hour (ex. 7 PM).
- Replace the battery when the LCD displays "LO".
- When the LCD screen displays „Err“ you could have exceeded the 150Kg limit.

## Calculate your BMI !

### (Body Mass Index)

BMI (Body Mass Index) is a number that is calculated using the weight and height of a person. It is a pretty accurate indicator of the fat in the body for most people. BMI does not measure the body fat in a direct way. BMI can be considered an alternative to measuring the body fat directly. What is more, calculating the BMI is a cheap and accessible way of monitoring the weight categories that could create health issues.

For women, an excellent BMI is between 21 and 23. Overweight is indicated by a BMI 20% bigger than the highest possible value (over 27.3), while obesity has values 40% bigger than the maximum ideal indicator (over 32.3).

For men, an excellent BMI is between 22 and 24. Overweight is indicated by a BMI 16% bigger than the highest possible value (over 27.8), while obesity has values 40% bigger than the maximum ideal indicator (over 31.1).

Calculating the BMI is recommended for persons over 20 years old. In case of children and adolescents calculating it is more complex and that is why different methods of calculating the BMI are applied.

### BMI INTERVALS - men

**Underweight:** < 20,7

**Normal weight:** 20,7 – 26,4

**Slightly overweight:** 26,4 – 27,8

**Overweight:** 27,8 – 31,1

**Obesity:** 31,1 – 45,4

**Extreme obesity:** 45,4 <

**MEN**

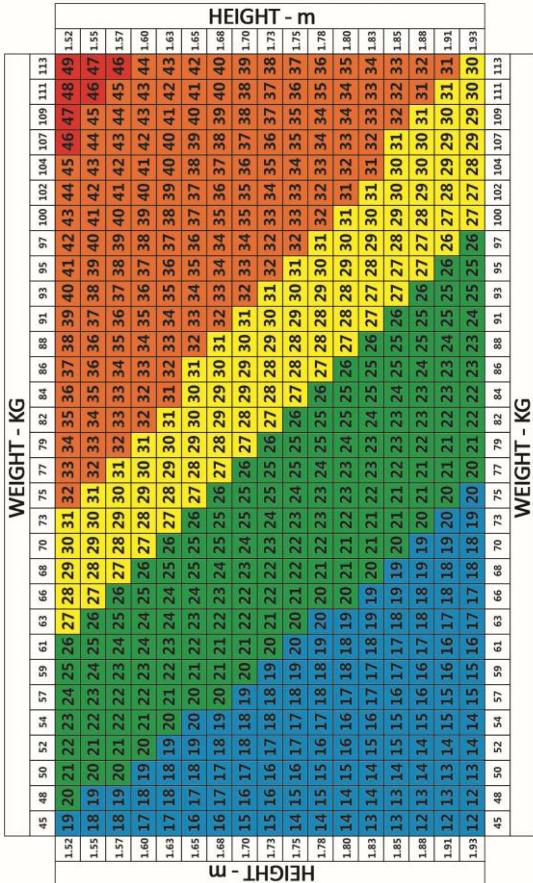
Table of optimal weight:

Women		Men	
Height (cm)	Weight (kg)	Height (cm)	Weight (kg)
148	43	155	60
150	45	160	63
153	47	163	65
155	49	165	66
160	52	168	69
163	53	170	71
165	55	173	73
168	57	175	74
170	60	178	77
173	62	180	78
175	64	183	81
178	67	185	83
180	70	188	86
183	72	190	88
185	73	195	90

CEP 1221 Specifications	
Measured weight	From 6Kg to 150Kg
Divisions	0,1 Kg
Display	LCD Screen
Measurement technology	High precision sensors
Measuring unit	Kg, lb, st.
Current supply	3V (1x CR 2032, battery included in pack)
CEP 1221 Characteristics	
Material	Secure glass
Device weight	1.40 Kg
Device dimensions	Diameter 330mm, Height 22mm
Package weight	1.60KG
Package dimensions	L340xW337xH30 mm
Color	Transparent
Storage temperature	From 0° C to 50° C
Temperature of use	From 10° C to 40° C

**CEP 1221**

**Recycling:**

This product is marked with the distinctive mark of selective sorting for electric and electronic waste. This symbol indicates the fact that this product must be accepted by a selective waste collection center, in the accord with the directive EU 2002/96/CE, so that it can be recycled or disassembled to reduce the impact to the environment.  
**Beware!** Electronic products, unsorted selectively, may affect the health of the environment and people because they may contain dangerous substances.

**BATTERIES:**

- Throwing the batteries in the garbage is forbidden. Take the used box to recycle centers.
- Incorrect use of batteries may cause leaks or corrosion.
- Do not mount the batteries incorrectly.
- Charging/hitting/opening/short-circuiting the batteries is dangerous.
- Do not let drain batteries inside the scale.
- Replace batteries when the scale is functioning improperly.
- Clean the liquid and replace the batteries if they have leaked inside the scale.
- The batteries inside the pack do not contain heavy metals.

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